




## CD-Imager

**The innovative robotic system for a professional image display on CD or DVD**

- Fully automatic compact system for creating DICOM patient CDs or DVDs
- Highly compatible with all digital DICOM modalities (multimodality)
- Individual labelling (practice/clinic logo)
- Easy integration of DICOM patient data
- Extremely cost effective due to quick printing times and low link consumption

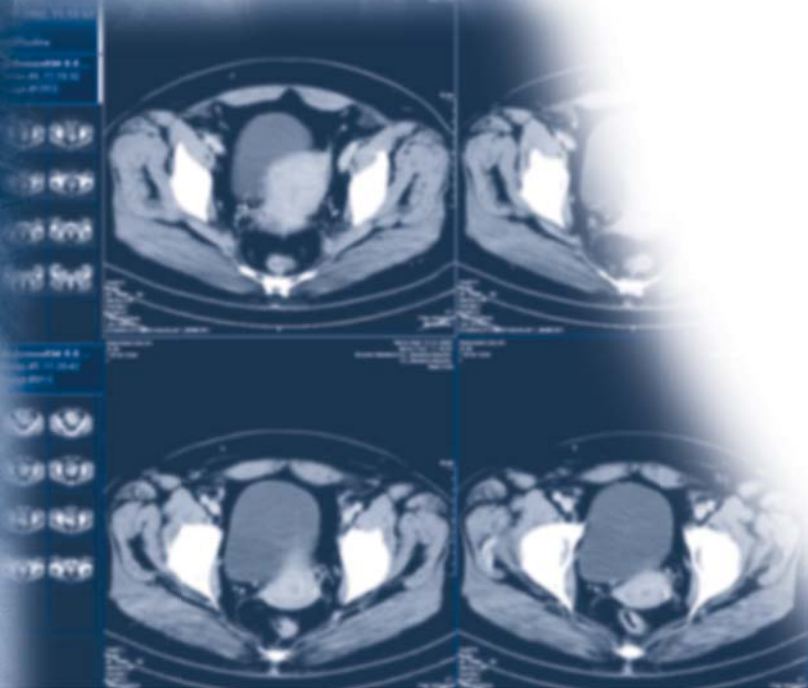

## DIGITAL IMAGES ARE OUR WORLD

The medigration CD Imager is a cost effective, fully automatic robotic system for creating CDs or DVDs in DICOM offline media format. It enables you to burn CDs or DVDs of image data and findings for referral physicians or patients in a way that is both

professional yet uncomplicated. The images can be displayed on any PC with the integrated viewer (ImageVision light) in original DICOM quality and can be imported to all PACS systems totally independent from the respective manufacturers.

### Extremely easy to use


With our medigration CD Imager it is possible to begin the burning process directly from the modalities or PACS stations within the practice network. As well as patient data, there is the option to transfer the submitted results from the RIS/HIS to HL7. With the easy to manage user interface, you have all running work processes in view all the time and can intervene when necessary. In the sense of corporate identity, our system enables individual labelling by integrating your practice/clinic logo into your print layout. We can incorporate two different logos on request.



Individual labelling

### Integrated DICOM viewer

As well as the standard functions (windowing, zoom, inverting, rotation, annotation) our integrated viewer ImageVision light offers the option of displaying automatic cut-lines in images with defined positional relations. In addition, two separate layouts can be set out in the comparison mode. After closing the DICOM viewer (ImageVision light), your desired website opens automatically (homepage of your clinic/practice).

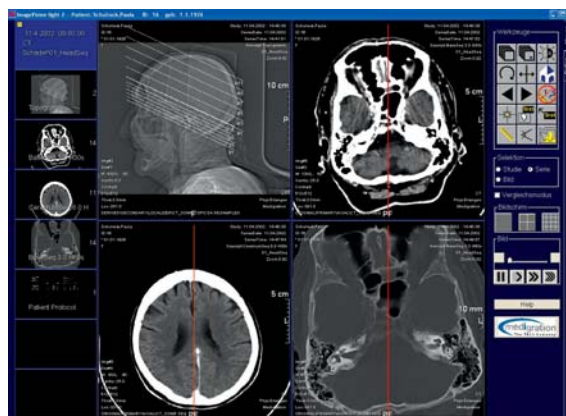

Ergonomic user interface ImageVision

### Technical Details

A high media presentation capacity (up to 300 discs), processing various media formats at once, an additional output compartment on the side and the "Perfect Print" function make the CD Imager Advanced a versatile and flexible production system. It offers professional and robust

printing technology, which is particularly suited for medium to high production output. In combination with Perfect Print, silk screen pre-printed media can also be subsequently printed.

The **money-saving** alternative to laser film!



Gesellschaft für  
medizinische Datenverarbeitung  
und Integration mbH  
Schuhstr. 30  
D-91052 Erlangen

Tel: + 49 (0) 91 31 - 6 90 87 - 40

Fax: + 49 (0) 91 31 - 6 90 87 - 50

eMail: info@medigration.de

web: www.medigration.de



For more information please visit  
www.medigration.de  
If you have any further questions  
please feel free to contact us!